

SHINING 3D  
DENTAL

# MetiSmile

3D Face Scanner

Beyond Smile

 **OralTech**  
DIGITAL




SHINING3DDENTAL.COM





# SAY HELLO TO MetiSmile

MetiSmile is the first face scanner developed and produced by SHINING 3D exclusively for dentistry. It can quickly capture the facial information to create 3D model and assist in clinical diagnosis with its advanced software.

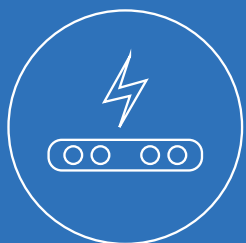


**Fast Scan Speed**



**Ortho Simulation**



**Mandibular Trajectory Tracking**

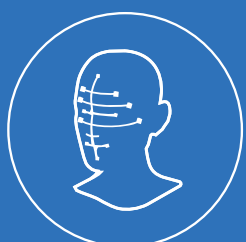

**Facial Features  
Measurement & Comparison**



**Automatic Alignment of  
Intraoral and Facial Data**




 **OralTech**  
DIGITAL







# POWERFUL HARDWARE

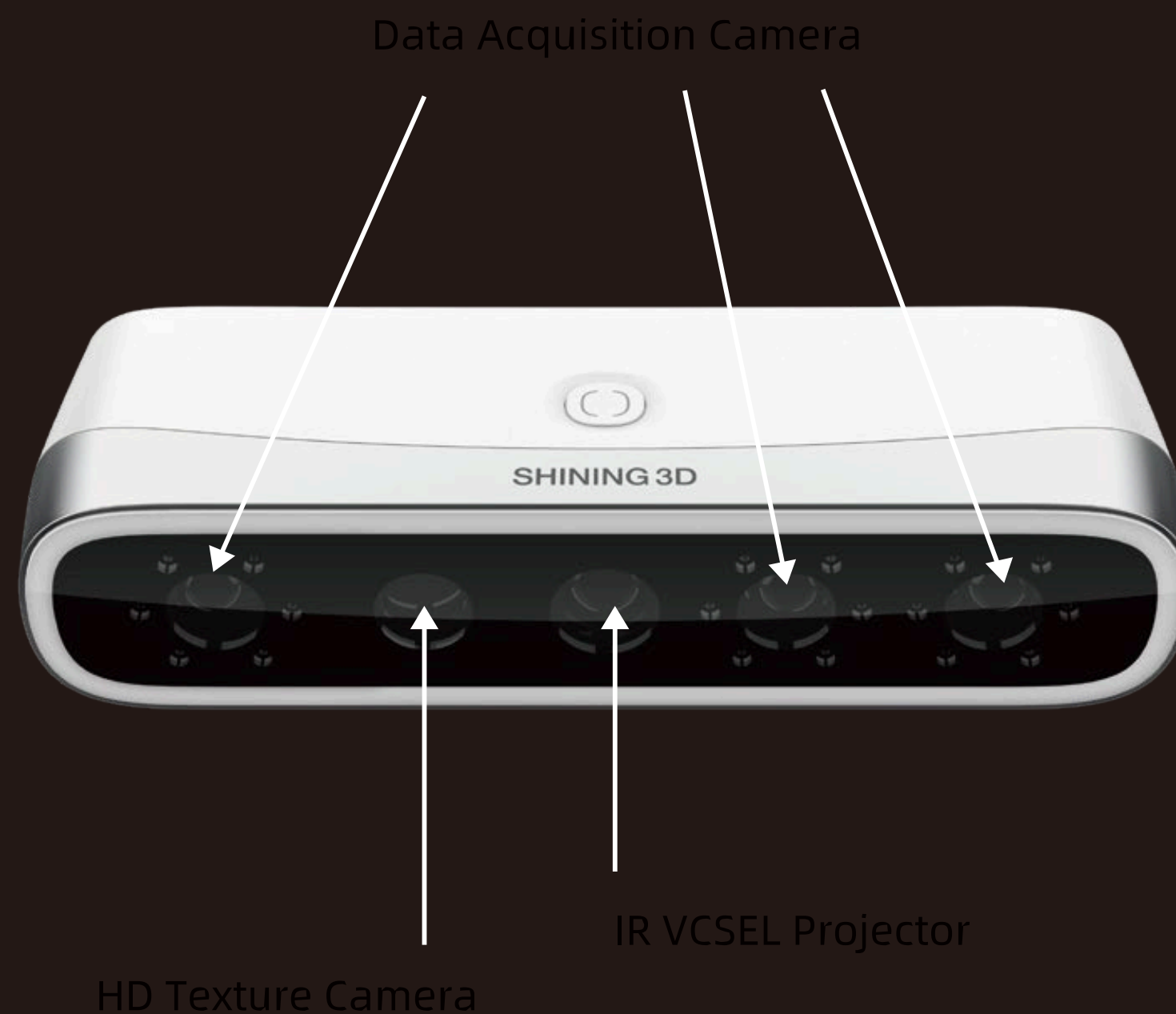

## Fast Scan Speed

In just 10 seconds this highly perceptive scanner can take photos of patients from multiple angles to simultaneously construct a 3D facial data.

## High Accuracy

Three 1.3 MP data acquisition cameras and one 5.0 MP HD texture camera produce scan accuracy within 50 $\mu$ m. MetiSmile also captures elevated details of the teeth.

## High-fidelity Texture

The exceptional texture camera of the MetiSmile can accurately record and display facial color that appears realistic to the patient.



# ADVANCED SOFTWARE


## Auto Alignment Between Facial Scan Data and Intraoral Scan Data

The software will align the facial and intraoral scans which displays the patients' facial and oral information in a comprehensive way. The scans are overlapped to provide a better view and perspective for digital smile design that assists in creating more aesthetic restorations for patients.



### Step 1

Use Intraoral scanner to capture the intraoral data.

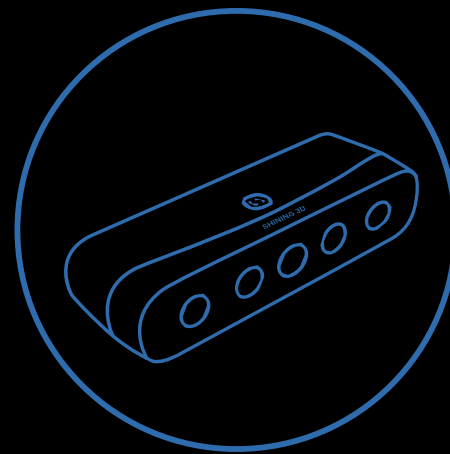

### Step 2

Use MetiSmile to capture the facial data.



### Step 3

Import the intraoral data to MetiSmile software for auto alignment.






Advanced Software

# ORTHO SIMULATION

Some of the features of the ortho simulation module include automatic lip extraction and one-click tooth segmentation, etc. The intelligent before and after orthodontic treatment display is enhanced by the integration of the facial scan information.


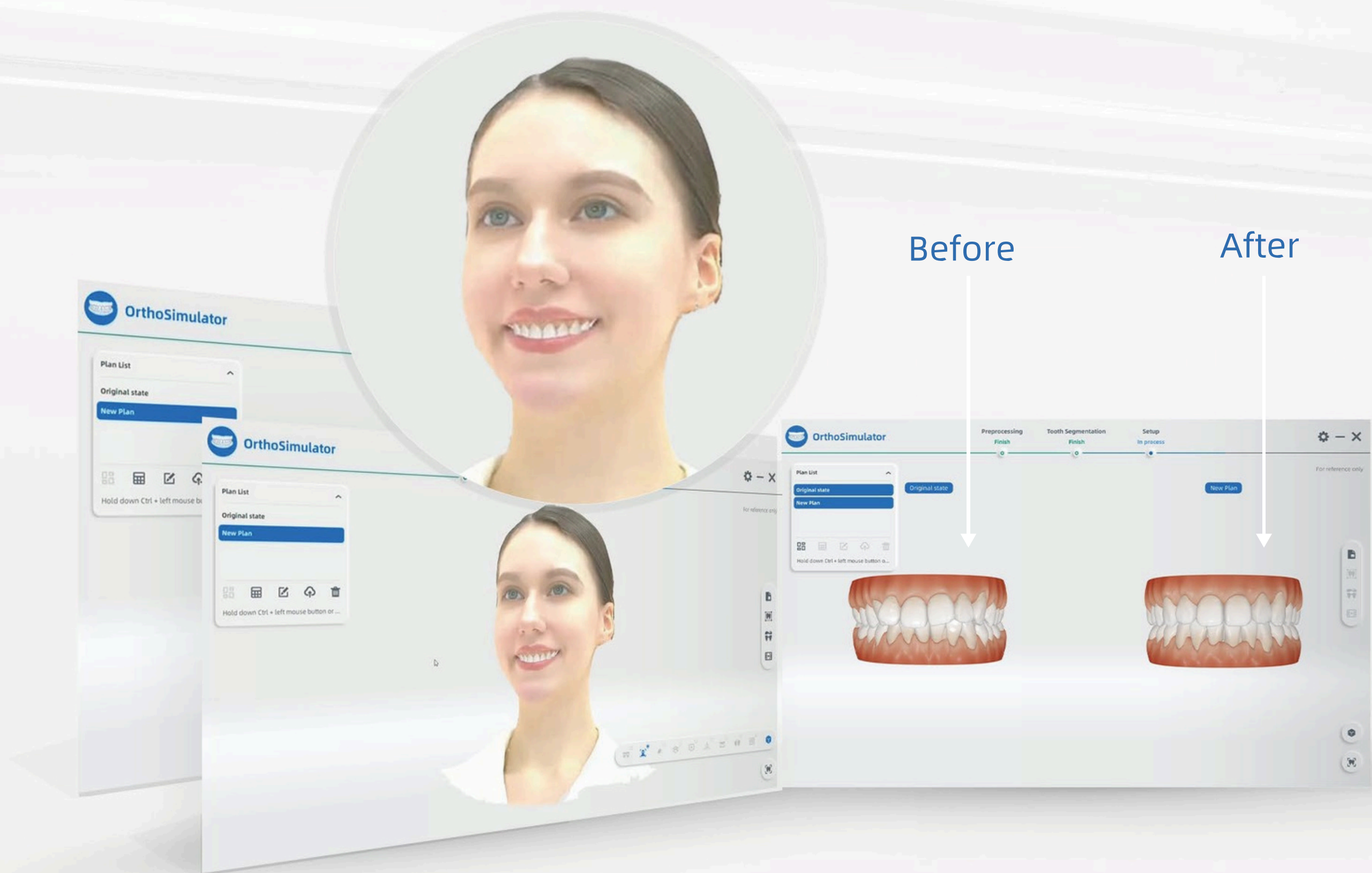




Advanced Software

# FACIAL FEATURE MEASUREMENT & COMPARISON

The intuitive software is capable of recognizing 3D facial features. It can measure the distance of facial data, and compare facial changes before and after treatment. These software tools help dentists diagnose and evaluate treatment more effectively. It also improves the patient's satisfaction with the treatment plan while adding visual presentation.

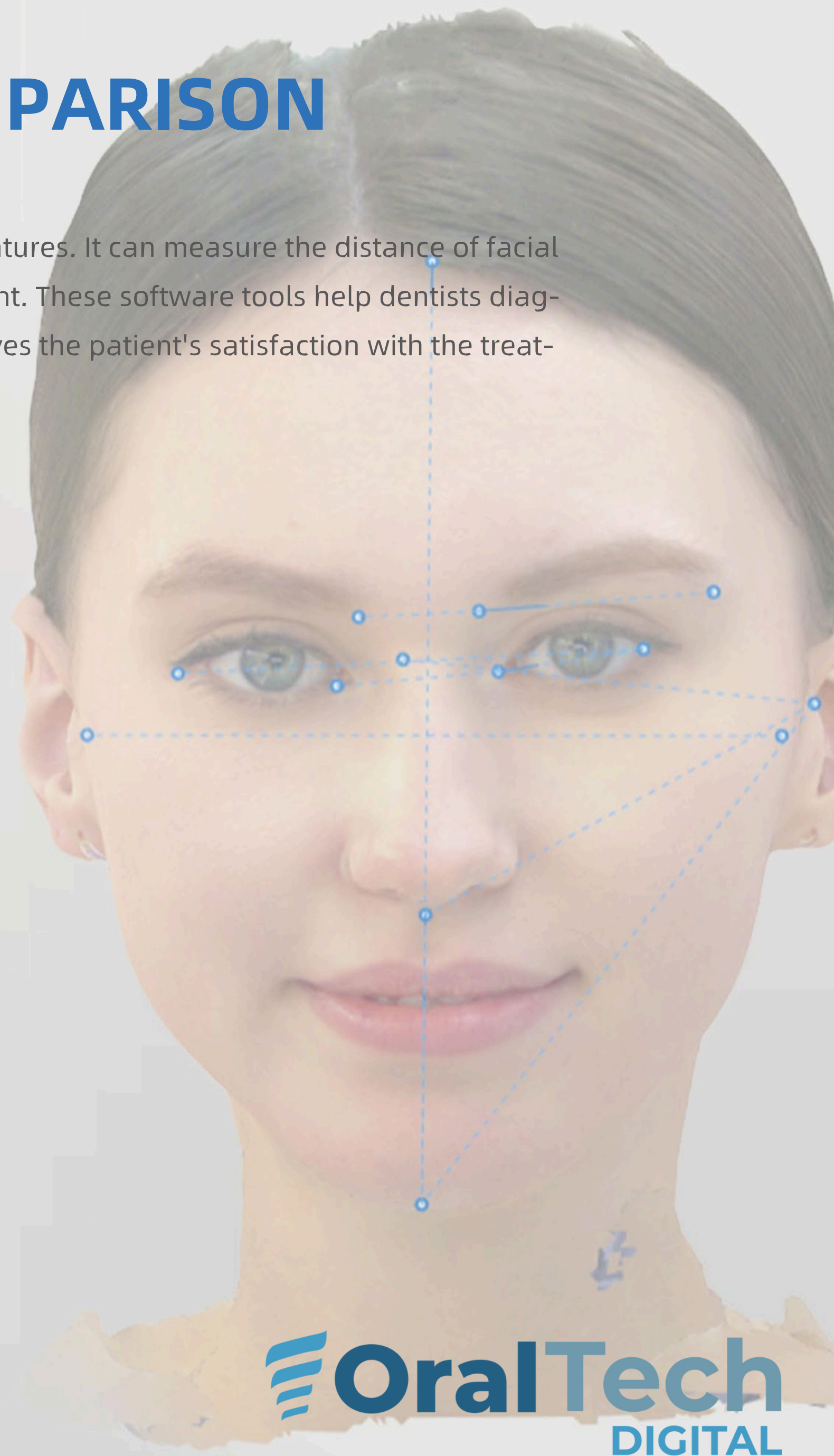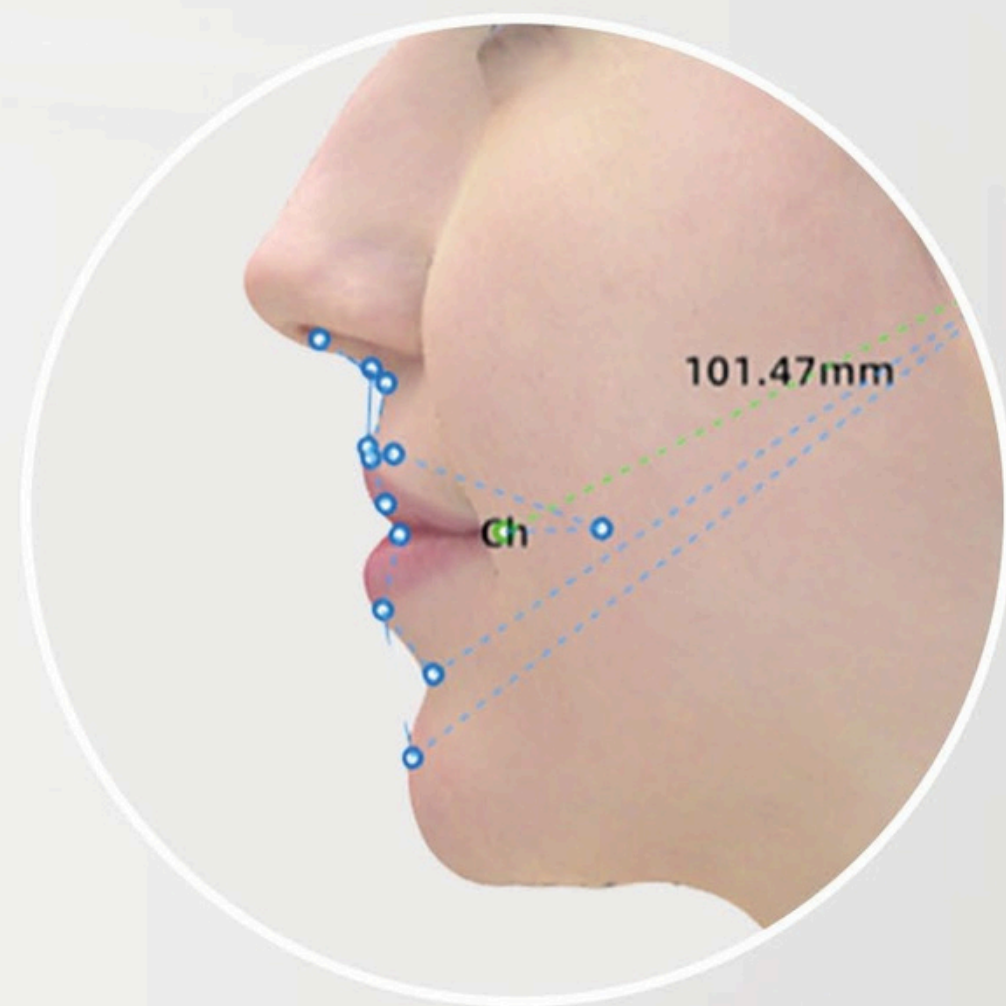



Advanced Software

# MANDIBULAR TRAJECTORY TRACKING

Using mandibular trajectory tracking, the dentist can acquire dynamic occlusion data. This module includes left and right lateral detection, centric, and open occlusion. It delivers detailed occlusion information for accurate diagnosis, design, and treatment.





# APPLICATIONS

The powerful and advanced software of the MetiSmile makes this product an indispensable tool for digital dental treatment, including maxillofacial surgery, implant and prosthesis, orthodontics, medical cosmetology, etc. It will take your clinic to the next level.



Maxillofacial surgery



Implant and prosthesis



Orthodontics



Medical cosmetology



Aesthetics restoration







# EXTRAORDINARY USER EXPERIENCE

- ❑ Eye-friendly flashless scan thanks to infrared technology.
- ❑ Automatic brightness adjustment ensures outstanding face texture.
- ❑ Guided operation throughout the scan workflow.
- ❑ Open system to export STL, OBJ and PLY.
- ❑ Elegant and compact, only 800 grams.

Handheld Mode

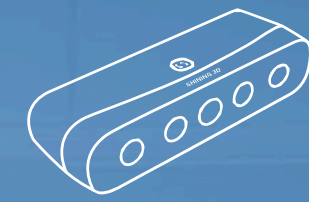

Fixed Mode

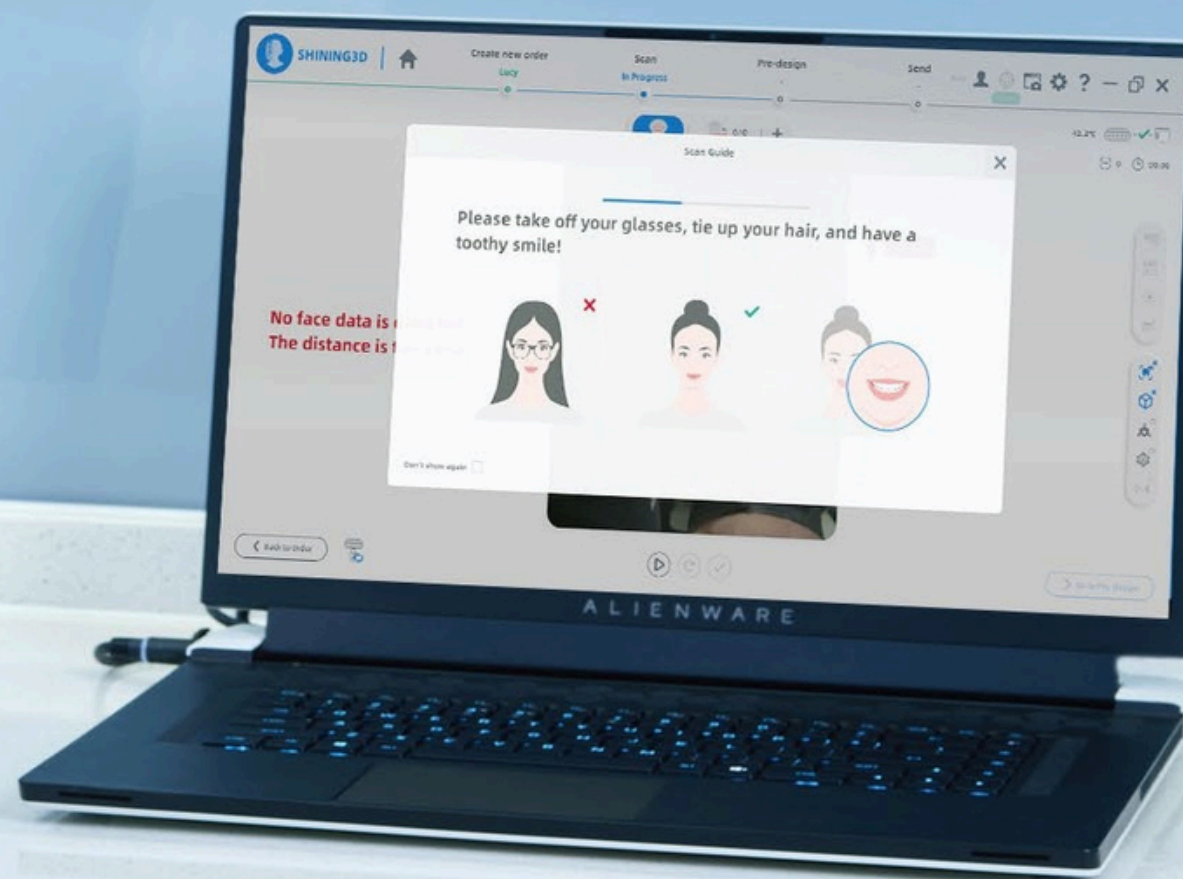

 **OralTech**  
DIGITAL



# Technical Specifications

## MetiSmile

Model	MetiSmile
Resolution	Data Acquisition Camera: 1.3 Mega Pixel HD Texture Camera: 5.0 Mega Pixel
Accuracy	50µm
Field of View	With working distance 500mm, the FOV is 210*270mm
Output Format	PLY, OBJ, STL
White LED color temperature	5500K
Dimension	215*50*75mm
Weight	800g
Power Supply	Input: AC100-220V~, 50/60HZ, 1.5A Output: DC12V, 7.0V

### Recommended PC Configuration


CPU	Intel Core i7-8700 or higher
Memory	16GB is the minimum, 32GB is highly recommended.
Hard Disk Drive	256GB SSD or above
Display Resolution	1920*1080, 60Hz or higher
Connector	USB 3.0
Graphic Card (GPU)	NVIDIA RTX 2060 6GB or higher
Operating System	Microsoft Windows 10 (64-bit) or later versions of Windows operating system



# ABOUT US

SHINING 3D provides fully integrated 3D digital dental solutions. Starting with obtaining scan data with 3D scanners for labs, intraoral scanners for clinics, to designing with dental CAD software. From case design, to 3D print, including: working models, orthodontics models, implant models, surgical guides, wax-ups, and partial frameworks.

SHINING 3D obtained CE, FDA international product standards, ISO13485 medical device quality management system certification, Authorized Economic Operator (AEO) certification; ISO/IEC 27001, 27701, 27017 and 27018 four IT security certifications.

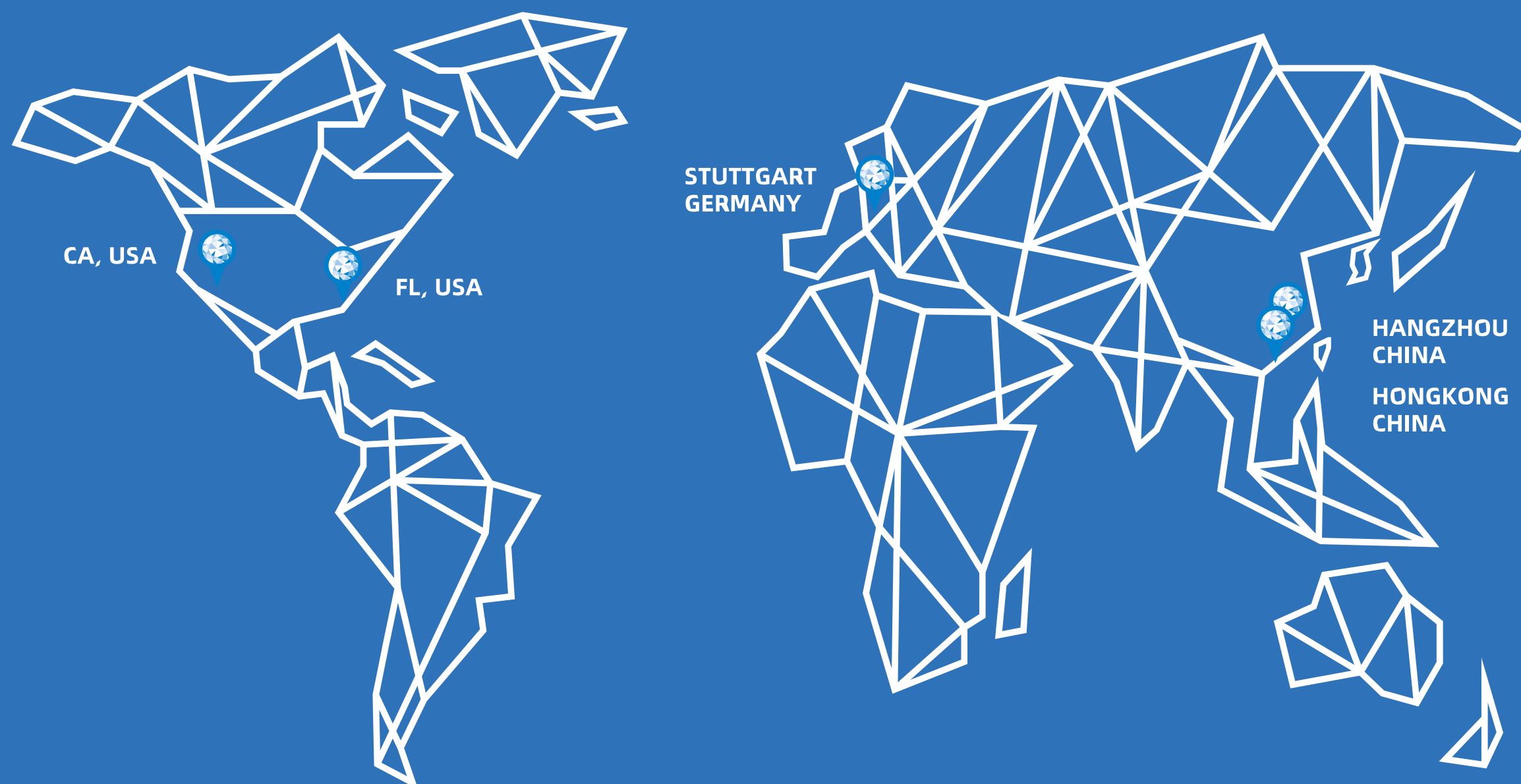

## Global Headquarters

SHINING 3D Tech. Co., Ltd.  
Address: No. 1398, Xiangbin Road, Wenyan,  
Xiaoshan, Hangzhou, Zhejiang, China, 311258  
Tel: +86 571 8299 9050

## EMEA Region

SHINING 3D Technology GmbH.  
Address: Breitwiesenstraße 28 70565  
Stuttgart, Germany  
Tel: +49-711 28444089

## Americas Region

ORALTECH LATAM  
AV. MANUEL J. CLOUTHIER 699  
ZAPOPAN JALISCO MEXICO  
TEL: 332561107